

BILL TO:

Account Name: _____

Account #: _____ P.O. #: _____

Contact: _____ Phone: _____

Address: _____

City: _____ State: _____ Zip: _____

SHIPPING INFORMATION: Same as "BILL TO"

Name: _____

Address: _____

City: _____ State: _____ Zip: _____

 UPS Red (overnight) UPS Blue (2 day air) UPS 3 Day UPS Ground Other _____**PRODUCT SELECTION: CCMi Mark III**Frame: DEFIANCE® III DEFIANCE®

Hinge:

 FourcePoint™ (ACL only) Female Fource™ (ACL only) Standard (D III) DropLock (CI only) Elbow (CI only) Adjustable OA

Strapping:

 ACL PCL CIOA Affected Compartment: Medial Lateral

Degrees of Correction: _____ (1° to 7°) (Standard= 3°)

FLAT COLORS:

- Matte Black Neon Pink Neon Yellow Neon Orange
 Black Red Purple Pink
 White Blue Light Blue Yellow
 Green Orange Tan Red, White, Blue

METALLIC COLORS:

- Black Red Purple Pink
 Navy Dark Blue Bright Blue Gold
 Green Orange White Silver

MULTIPLE COLORS:

Thigh Frame: _____

Calf Frame: _____

Fade (flat colors only): Yes No**GRAPHICS (extra charge):**

- Tiger Red Tiger Orange Carbon Fiber
 Snake Camouflage
 Stars & Strips (flat) Stars & Strips (metallic)

DONJOY

#1 in Knee Bracing

CCMI Mark III1430 Decision Street
Vista, CA 92081-8553
USATel: (800) 336-6569
www.djglobal.com**Fax To:** (800) 936-6569**Email To:** customercare@djglobal.com**FOR OFFICE USE ONLY:**

Order #: _____

Brace #: _____


PATIENT INFORMATION:

Name: _____

Age: _____ Height: _____ Weight: _____

Knee Measurement: Right Left Reverse

Measured By: _____

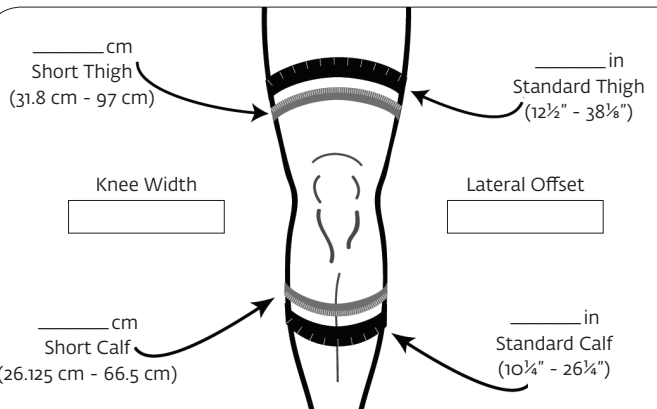 New Brace Remeasurement/Repair Refurbish

Note: Check both short and standard measurements for ski boot option. Best to measure with Ski Boot On.

FEATURES (*standard option on DEFIANCE® III):

- Suede Liners/Pads
 Chamois Liners/Pads
 Swooping Thigh (not available on short thigh)*
 Swiveling Strap Tabs*
 Metal Caps - PowerCaps (recommended for contact sports)
 Reinforced Carbon Composite Frame
 Legend Strap Ends
 Installed Extension Stop: 0° 10° 20° 30° 40° (circle one)
 Installed Flexion Stop: 45° 60° 75° 90° (circle one)
 Neoprene Suspension Strap Kit

ACCESSORIES (extra charge):

- High Strength Stainless Steel Gears Air Condyle Kit
 Full Extension Lockout Stop Foam Condyle Pads
 Lycra® Suspension Sleeve Armor Impact Guard
 Neoprene Suspension Sleeve Pinch Guard (Calf)
 Sports Brace Cover Atrophy Pad Kit
 Sili-Grip Strap Pad Kit Tru-Pull® Advanced Attachment


For further instructions, please refer to the CCMII Mark III instruction video and instruction card included with the instrument.

CAUTION!

FEDERAL LAW (U.S.A.) RESTRICTS THIS DEVICE TO SALE BY THE ORDER OF A LICENSED HEALTH CARE PROFESSIONAL

Custom Contour Measuring Instrument - Mark III Use Instructions

NOTICE: For optimum long term fit of the DEFIANCE® or DEFIANCE® III, measurements should be taken when there is no edema or atrophy present.

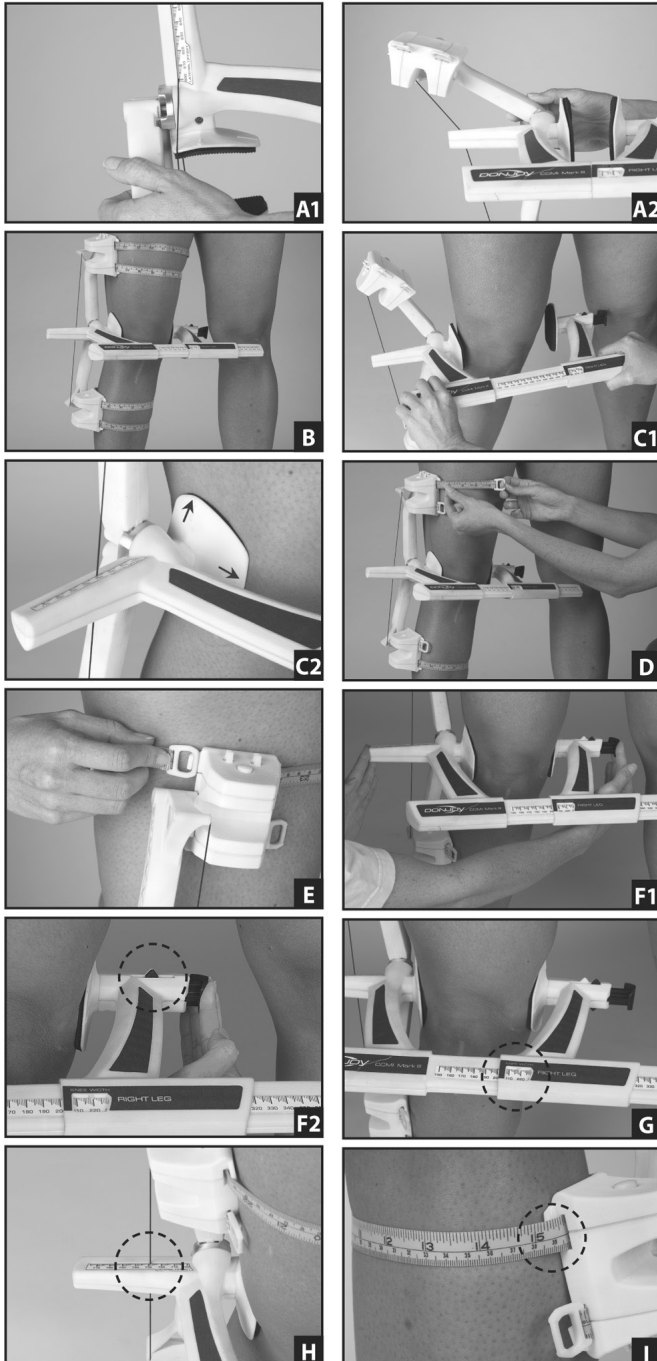

APPLICATION INFORMATION:

1. To unstow, distract and rotate swing arms as shown in figures (A1 and A2).
2. Have the patient stand feet shoulder width apart and knees flexed slightly. Do not allow the knees to flex beyond 15 degrees. Apply the instrument directly to the bare leg. Note: Taking measurements from the uninvolved leg is not recommended because of differences in leg dimensions due to dominance of one leg over the other. See figure B.
3. Distract knee width caliber and position on the patient's knee. Arrows on condyle shells should align with the top of the knee cap and slightly posterior to midline. See figures (C1 and C2).
4. For standard thigh and calf, extend tape measure marked in Inches located distal to knee. Wrap around leg and secure to hook on tape module. Press tape release button to pull tape snug to leg; ensure it is parallel to floor. For short thigh or calf, follow same procedure using Metric measurements. See figure (D and E).
5. To establish proper knee compression for knee width measurement, support lateral side of instrument and press load plunger on medial condyle until triangles are aligned, as shown in figures (F1 and F2).
6. Record knee width on order form. See figure (G.)
7. Record lateral offset on order form. See figure (H).
8. Record thigh and calf circumference on order form. See figure (I).